



Procedure to access the power supply and TEC controllers on the Marine Aeri

1. Open the back door
2. Remove the isolation cover (2 screws), this is thin aluminum cover with a isolation foam inside.

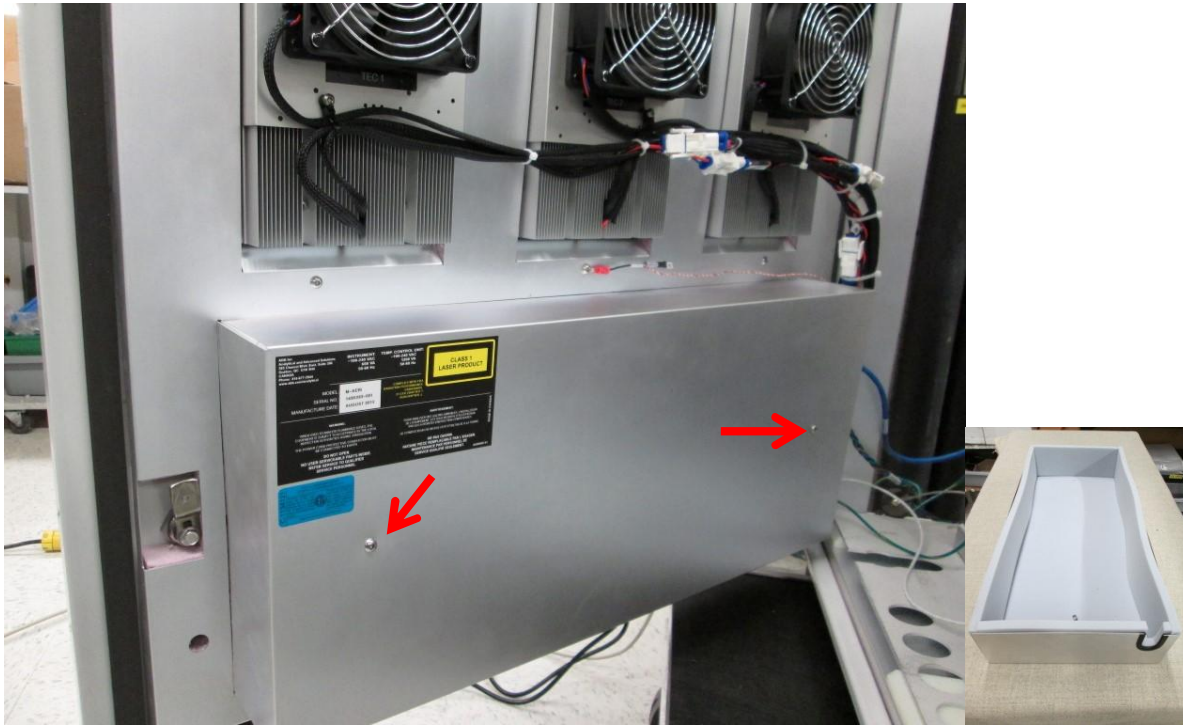

3. Remove the power supply and TEC controller cover (12 screws), this stainless steel cover has 12 screw that are mounted with an O-ring to facilitate the replacement of the panel. Pay attention to the cables on the right side when taking it out.

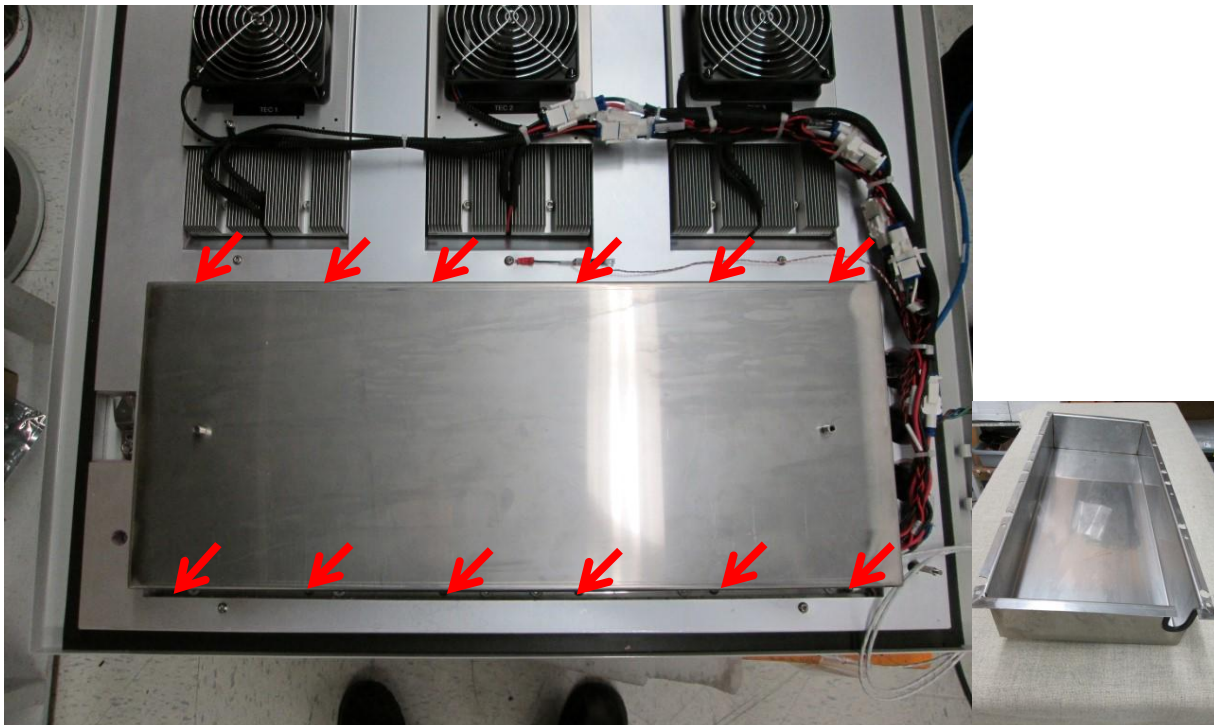

4. Once the Power supply and TEC cover removed, you will have access to the power supply and the TEC controller. This view show the door removed from the main frame. Don't take off the door.

